

AvaSpec-HERO SensLine

The AvaSpec-HERO is the top of the line spectrometer!
Based on our High Sensitivity Compact (HSC) optical bench (f=100mm; NA=0.13) and a 1024x58 backthinned CCD detector, it offers the best of both worlds: high sensitivity and resolution!

The instrument is equipped with thermo-electric cooling, enabling long integration times in low light applications. In conjunction with our AS7010 electronics, including a high-end AD converter, noise is kept to a minimum, which offers you an excellent Signal to Noise and Dynamic Range performance.

A selection of gratings and slits offers you the flexibility of configuring the instrument for a wide range of applications in the 200-1160 nm range.

From low light fluorescence applications to demanding Raman applications, the AvaSpec-HERO is your ideal companion.

With the high-speed USB3.0 and Gigabit Ethernet communication interface, the connection to your computer is fast and simple.

Of course the digital IO ports enabling external triggering, control of shutters, and pulsed light sources from the Avantes line of instruments are available as well.

The Avaspec-HERO is standard equipped for use with replaceable slits, offering optimal flexibility for a variety of applications. The combination of all the above makes the AvaSpec-HERO your ideal companion for all your spectroscopic measurements.

AvaSpec-HERO



Technical Data

Optical Bench	HSC Symmetrical Czerny-Turner, 100 mm focal length, NA: 0.13
Wavelength range	200-1160 nm
Resolution	0.2-7 nm, depending on configuration (see table)
Stray-light	0.5%, depending on the grating
Sensitivity	445,000 counts/ μ W per ms integration time
Detector	CCD array image sensor with one stage TE Cooled, 1024 pixels
Temperature cooled CCD	Max. $\Delta T = 30$ °C versus ambient
Signal/Noise	1200:1
Dynamic Range	40.000
AD converter	16-bit, 250 kHz
Integration time	5.2 ms- 60 sec
Interface	USB 3.0 high-speed, 5 Gbps Gigabit Ethernet 1 Gbps
Digital IO	HD-26 connector, 2 Analog in, 2 Analog out, 3 Digital bidirectional, trigger, sync., strobe, laser.
Sample speed with on-board averaging	5.2 ms/scan
Data transfer speed	5.2 ms/scan (USB3 and ETH)
Power supply	12VDC, 1.5A
Dimensions, weight	185 x 161 x 185mm, 3500 grams

Distributor



nordics@amstechnologies.com
azpect.amstechnologies.com

Contact us 

The new AvaSpec-HERO is the answer for those who are in need of high resolution and high sensitivity!

Grating Selection Table for AvaSpec-HSC1024x58TEC-EVO

Use	Useable range (nm)	Spectral range (nm)	Lines/mm	Blaze (nm)	Order code
UV/VIS/NIR	200-1160	770-760*	300	300	HSC0300-0.30
UV/VIS/NIR	250-1160	770-760*	300	420	HSC0300-0.42
VIS/NIR	250-1160	577-553	400	550	HSC0400-0.55
UV/VIS	250-850	373-340*	600	400	HSC0600-0.40
VIS/NIR	250-1160	373-340*	600	650	HSC0600-0.65
VIS/NIR	500-1160	268-220*	830	900	HSC0830-0.90
UV/VIS	200-1160	182-130*	1200	400	HSC1200-0.40
VIS/NIR	500-1050	182-130*	1200	750	HSC1200-0.75
UV/VIS	200-580	84-61*	2400	270	HSC2400-0.27

* depends on the starting wavelength of the grating; the higher the wavelength, the bigger the dispersion and the smaller the range to select.

Resolution Table (FWHM in nm) for AvaSpec-HSC1024x58TEC-EVO

Grating (lines/mm)	Slit size (μm)					
	10	25	50	100	200	500
300	1.70	1.80	2.30	3.40	6.50	14.0
400	1.30	1.45	1.60	2.60	5.10	12.0
600	0.75	0.85	1.10	1.70	3.40	7.50
830	0.50	0.60	0.70	1.25	2.30	5.00
1200	0.32	0.40	0.48	0.80	1.45	3.50
2400	0.17	0.30	0.36	0.50	0.80	1.75

* Above values are average values. Due to optical properties resolution will be better in the lower wavelengths than in the higher wavelength range.

Ordering Information

AvaSpec-HSC1024x58TEC-EVO

- AvaSpec-HERO; High sensitivity fiber optic spectrometer, HSC 100mm bench design, 1024x58 pixel back illum TE cooled CCD detector, high-speed USB 3.0 and ETH interface, including AvaSoft-Basic, USB interface cable, specify grating, wavelength range and options

Options

- SLIT-XX-RS** • Replaceable slit with SMA connector, specify slit size XX=10, 25, 50, 100, 200 or 500 μm.
- SLIT-XX-RS-FCPC** • As SLIT-XX-RS, but with FC/PC connector
- SLITKIT-SMA** • Slit kit containing 25, 50, 100, 200 or 500 μm slits, and the tools to replace the slit. SMA-connectors
- SLITKIT-FCPC** • As SLITKIT-SMA, but with FC/PC connectors
- OSF-YYY-3** • Order sorting filter for reduction of 2nd order effects, 3 mm thick, please specify YYY=305, 395, 475, 515, 550, 600 nm
- OSC-HSC300** • Order sorting coating for use with grating HSC0300-xx
- OSC-HSC600** • Order sorting coating for use with grating HSC0600-xx and HSC0400-xx



nordics@amstechnologies.com
azpect.amstechnologies.com

Contact us 

AvaSpec-HS2048XL-EVO SensLine High UV and NIR Sensitivity Back-thinned CCD Spectrometer

For high sensitivity applications where high resolution is not of paramount concern, the AvaSpec-HS2048XL-EVO is an exceptional instrument. Featuring Avantes' HS optical bench which has a full 0.22 numerical aperture for superior throughput, the AvaSpec-HS2048XL has a back-thinned CCD detector with 2048 pixels measuring 14X500 microns.

Unlike many back-thinned CCD spectrometers, which have two dimensional arrays the HS2048XL has large monolithic pixels with exceptional efficiency in the UV, from 200-400 nm, and the NIR, from 950-1160 nm, while retaining sensitivity in the visible range. The unique optical design features torroid collimating and focusing mirrors to control image magnification and enhance efficiency. The instrument also features an electronic shutter, which enables integration times as low as 2 microseconds. For configurations, which require second order filtering, order-sorting filters are avail-

able. The AvaSpec-HS2048XL is available with a wide range of slit sizes, gratings and may be configured with SMA or FC/ PC fiber-optic entrance connectors.

The AvaSpec-HS2048XL-EVO uses the AS7010 electronics board offering USB3 (10 times faster than USB2), Gigabit Ethernet and better signal processing.

AvaSpec-HS2048XL-EVO



Technical Data

Optical Bench	High-sensitivity asymmetrical design, 37.5 mm focal length; NA - 0.22, f/2.27
Wavelength range	200 - 1160 nm
Resolution	1 - 20 nm, depending on configuration (see table)
Stray-light	< 1 %
Sensitivity	1,250,000 counts/ μ W per ms int. time
UV Quantum efficiency	60% (200-300 nm)
Detector	Back-thinned CCD image sensor 2048 pixels
Signal/Noise	525:1
AD converter	16-bit, 1 MHz
Integration time	2 μ s - 600 seconds
Interface	USB 3.0 high-speed, 5 Gbps Gigabit Ethernet, 1 Gbps
Sample speed with on-board averaging	2.44 ms /scan
Dynamic Range	14.900
Data transfer speed	2.44 ms /scan (USB3)
Digital IO	HD-26 connector, 2 Analog in, 2 Analog out, 3 Digital in, 12 Digital out, trigger, synchronization
Power supply	Default USB power, 700 mA. or external 12VDC, 360 mA
Dimensions, weight	175 x 165 x 85 mm, 1,950 kg

Grating Selection Table for AvaSpec-HS2048XL-EVO

Use	Useable range (nm)	Spectral range (nm)	Lines/mm	Blaze (nm)	Order code
UV/VIS/NIR	200-1160	900	500	330	HS500-0.33
UV/VIS	200-660	440	1000	250	HS1000-0.25
UV	200-850	520	600	300	HS600-0.30
UV/VIS	200-850	520	600	400	HS600-0.40
UV/VIS	300-1160	860	500	560	HS500-0.56
VIS	360-1000	500	600	500	HS600-0.50
NIR	500-1050	500	600	750	HS600-0.75
VIS	350-850	460	900	550	HS900-0.55
VIS	400-722	322	1200	500	HS1200-0.5
NIR	600-1160	500	600	1000	HS600-1.0
NIR	600-1160	350	830	900	HS830-0.9
NIR	750-990	240	1200	1000	HS1200-1.0

Resolution Table (FWHM in nm) for AvaSpec-HS2048XL-EVO

Grating (lines/mm)	Slit size (μm)					
	10	25	50	100	200	500
500	2.6	4.5	5.5	6.5	10.0	22.0
600	2.2	3.8	4.5	5.5	7.5	18.0
830*	2.1	3.6	4.0	5.0	7.0	15.0
900*	2.0	3.5	3.8	4.8	6.8	14.5
1000*	1.9	3.3	3.6	4.6	6.6	14.0
1200*	1.8	3.0	3.3	4.3	6.2	13.5

* theoretical values

Ordering Information

AvaSpec-HS2048XL-EVO

- High-sensitivity fiber-optic Spectrometer, 2048 large 500 μm pixel back-thinned CCD detector, USB powered, high-speed USB3.0 and ETH interface, incl. AvaSoft-Basic, USB interface cable. Specify grating, wavelength range and options

PS-12V/1.0A

- External power supply, needed for use in ETH mode

Options

SLIT-XX	• Slit size, please specify XX = 10, 25, 50, 100, 200 or 500 μm
OSF-YYY	• Order-sorting filter for reduction of 2nd order effects, 1 mm thick, please specify YYY= 305, 385, 475, 515, 550 or 600 nm
OSC-HS500	• Order-sorting coating with 350 and 600 nm long-pass filter for HS500 gratings in AvaSpec-HS
OSC-HS600	• Order-sorting coating with 350 and 600 nm long-pass filter for HS600 gratings in AvaSpec-HS
OSC-HS900	• Order-sorting coating with 600 nm long-pass filter for HS900 gratings in AvaSpec-HS
OSC-HS1000	• Order-sorting coating with 350 nm long-pass filter for HS1000 gratings in AvaSpec-HS
FCPC	• FC/PC fiber optic connector

The AvaSpec-HS2048XL-EVO is ideally suited for diffuse reflection measurements (UV, VIS, NIR) and fluorescence.

Contact us

nordics@amstechnologies.com
azpect.amstechnologies.com

Distributor
AZPECT
part of
amstechnologies
where technologies meet solutions

AvaSpec-ULS2048LTEC SensLine Thermo-Electric Cooled Fiber-optic Spectrometer

Long integration times in general are equivalent to higher dark noise. Avantes Thermo-Electric Cooled (TEC) spectro-meters systems overcome this problem


by cooling the detector. These instruments are equipped with triple stage cooling, keeping your detector at optimal 5 degrees Celsius (maximum -35°C difference from ambient temperature).

The detector cooling provides a significantly lower and more stable dark baseline and PRNU level. Dark noise is reduced by a factor of 2-3. This allows the ULS2048LTEC to be used in very low light conditions, such

as fluorescence and Raman measurements. If needed, integration times of more than 5 seconds are possible.

The AvaSpec-ULS2048LTEC has an integrated temperature regulator, USB2.0 high-speed interface and two cooling fans to actively ventilate the heat sink of the Peltier cooling elements. The spectrometer power supply is integrated into the housing.

AvaSpec-ULS2048LTEC

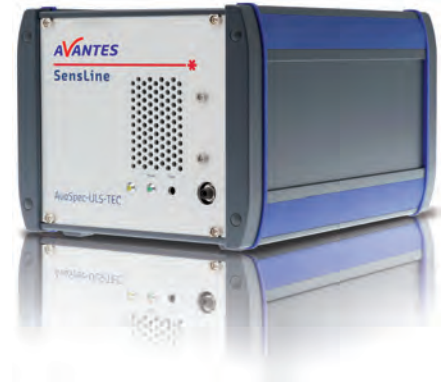

Technical Data

Optical Bench	ULS Symmetrical Czerny-Turner, 75 mm focal length
Wavelength range	200-1100 nm
Resolution	0.06 -20 nm, depending on configuration (see table)
Stray-light	0.04-0.1%, depending on the grating
Sensitivity	470,000 counts/μW per ms integration time
Detector	CCD linear array, 2048 pixels
Temperature cooled CCD	Max. ΔT = -30°C versus ambient (optimal setting: 5°C)
Time to stabilize	4 minutes
Dark baseline improvement @ ΔT=-35°C and it>5 sec	> Factor 6
PRNU improvement @ ΔT=-35°C and it>5 sec	> Factor 8
3-stage Peltier cooling internal Power supply @ ΔT=-35°C	5VDC, 3.0A
Signal/Noise	300:1
AD converter	16-bit, 2 MHz
Integration time	1.11 ms - 10 minutes
Interface	USB 2.0 high-speed, 480 Mbps RS-232, 115.200 bps
Sample speed with store to RAM	1.1 ms /scan
Data transfer speed	1.8 ms /scan (USB2) 430 ms/scan (RS-232)
Digital IO	HD-26 connector, 2 Analog in, 2 Analog out, 3 Digital in, 12 Digital out, trigger, sync.
Power supply	100-240 VAC, 50W
Dimensions, weight	250 x 179 x 144 mm, 3.6 kg

Distributor



nordics@amstechnologies.com
azpect.amstechnologies.com

Contact us

Our TEC-spectrometers are kept at a steady 5°C for maximum precision

Grating Selection Table for AvaSpec-ULS2048LTEC

Use	Useable range (nm)	Spectral range (nm)	Lines/mm	Blaze (nm)	Order code
UV/VIS/NIR	200-1100**	900**	300	300	UA
UV/VIS/NIR	200-1100**	900**	300	300/1000	UNA-DB
UV/VIS	200-850	520	600	300	UB
UV	200-750	250-220*	1200	250	UC
UV	200-650	165-145*	1800	UV	UD
UV	200-580	115-70*	2400	UV	UE
UV	200-400	70-45*	3600	UV	UF
UV/VIS	250-850	520	600	400	BB
VIS/NIR	300-1100**	800**	300	500	VA
VIS	360-1000	500	600	500	VB
VIS	300-800	250-200*	1200	500	VC
VIS	350-750	145-90*	1800	500	VD
VIS	350-640	75-50*	2400	VIS	VE
NIR	500-1050	500	600	750	NB
NIR	500-1050	220-150*	1200	750	NC
NIR	600-1100	350-300	830	800	SI
NIR	600-1100**	500**	300	1000	IA
NIR	600-1100	500	600	1000	IB

* depends on the starting wavelength of the grating; the higher the wavelength, the bigger the dispersion and the smaller the range to select.

** please note that not all 2048 pixels will be used for the useable range

Resolution Table (FWHM in nm) for AvaSpec-ULS2048LTEC

Grating (lines/mm)	Slit size (µm)					
	10	25	50	100	200	500
300	1.0	1.4	2.5	4.8	9.2	21.3
600	0.40-0.53*	0.7	1.2	2.4	4.6	10.8
830	0.32	0.48	0.93	1.7	3.4	8.5
1200	0.20-0.28*	0.27-0.38*	0.52-0.66*	1.1	2.3	5.4
1800	0.10-0.18*	0.20-0.29*	0.34-0.42*	0.8	1.6	3.6
2400	0.09-0.13*	0.13-0.17*	0.26-0.34*	0.44-0.64*	1.1	2.7
3600	0.06-0.08*	0.10	0.19	0.4	0.8	1.8

* depends on the starting wavelength of the grating; the higher the wavelength, the bigger the dispersion and the better the resolution

Ordering Information

- Thermo-Electric Cooled Fiber-optic Spectrometer, 75 mm Ultra-Low Stray-light AvaBench, 2048L pixel 3-stage TE-cooled and regulated CCD detector, USB2 high-speed interface, incl. AvaSoft-Basic, USB cable, desktop housing.
- Specify grating, wavelength range and options

Options

DUV	• Deep-UV detector coating >150 nm
DCL-UV/VIS-200	• Detector Collection Lens to enhance sensitivity, Quartz, 200-1100 nm
SLIT-XX	• Slit size, please specify XX = 10, 25, 50, 100, 200 or 500 µm
OSF-YYY	• Order-sorting filter for reduction of 2nd order effects, please specify YYY= 305, 395, 475, 515, 550 or 600 nm
OSC	• Order-sorting coating with 600 nm long-pass filter for BB (>350 nm) and VB gratings, recommended with OSF-305
OSC-UA	• Order-sorting coating with 350 and 600 nm long-pass filter for UA, VA gratings
OSC-UB	• Order-sorting coating with 350 and 600 nm long-pass filter for UB or BB (<350 nm) gratings
-FCPC	• FC/PC fiber optic connector
-RS	• Replaceable slit (recommended)

AvaSpec-ULS2048LTEC-USB2

nordics@amstechnologies.com
azpect.amstechnologies.com

Contact us

Distributor

AZPECT
part of
AMSTECHNOLOGIES

where technologies meet solutions

NEW

AvaSpec-ULS2048x64-EVO SensLine High UV and NIR Sensitivity Spectrometer

Alongside the cooled AvaSpec-ULS2048x64TEC-EVO with low-noise detector, Avantes also offers the more cost-effective, uncooled AvaSpec-ULS2048x64-EVO. With its standard 2048x64 backthinned detector, this spectrometer is perfect for less demanding applications in the UV and NIR range.

For applications that require integration times lower than 2 seconds, the cooling option is often not needed. For example, this uncooled AvaSpec-ULS2048x64-EVO has an established track record in various DOAS applications all over the world because of its high UV response and 0.9 mm detector height that enables detecting the wavelengths of interest.

Options include an order-sorting filter, to reduce second-order effects and purge ports for deep-UV measurements. The AvaSpec-ULS2048x64-EVO comes with a wide range of slit sizes, gratings and can be configured with SMA or FC/PC fiber-optic entrance connectors.

The AvaSpec-ULS2048x64-EVO uses the AS7010 electronics board offering USB3 (10 times faster than USB2), Gigabit Ethernet and better signal processing.

Connection to your PC is handled via USB3-connection or Ethernet, delivering a scan every 2 milliseconds. The instrument comes complete with AvaSoft-basic software, USB cable and an extensive manual.

AvaSpec-ULS2048x64-EVO



Technical Data

Optical bench	ULS, Symmetrical Czerny-Turner, 75 mm focal length
Wavelength range	200-1160 nm
Resolution	0.09-20 nm, depending on configuration (see table)
Stray light	< 1%, depending on the grating
Sensitivity	650,000 counts/μW per ms int. time
Detector	Back-thinned CCD image sensor 2048x64 pixels (height: 0.89 mm)
Signal/noise	450:1
AD converter	16-bit, 1.33 MHz
Integration time	2.4 ms-25 seconds
Interface	USB 3.0 high-speed, 5 Gbps Gigabit Ethernet 1 Gbps
Sample speed with on-board averaging	2.4 ms/scan
Readout noise	7.5 cnt RMS
Dark noise	11.5 cnt RMS
Dynamic range	6100
Data transfer speed	2.4 ms/scan (USB3)
Digital IO	HD-26 connector, 2 Analog in, 2 Analog out, 13 Digital bidirectional, trigger, sync., strobe, laser
Power supply	Default USB power, 885 mA. Or external 12VDC, 420 mA
Dimensions, weight	177 x 127 x 44,5 mm (1 channel), 1180 grams

Distributor

AZPECT
part of
amSTECHNOLOGIES
where technologies meet solutions

nordics@amstechnologies.com
azpect.amstechnologies.com

Contact us

Grating Selection Table for AvaSpec-ULS2048x64-EVO

Use	Useable range (nm)	Spectral range (nm)	Lines/mm	Blaze (nm)	Order code
UV/VIS/NIR	200-1160**	960**	300	300	UA
UV/VIS/NIR	200-1100**	900**	300	300/1000	UNA-DB
UV/VIS	200-850	520	600	300	UB
UV	200-750	250-220*	1200	250	UC
UV	200-650	165-145*	1800	UV	UD
UV	200-580	115-70*	2400	UV	UE
UV	200-400	70-45*	3600	UV	UF
UV/VIS	250-850	520	600	400	BB
VIS/NIR	300-1160**	860**	300	500	VA
VIS	360-1000	500	600	500	VB
VIS	300-800	250-200*	1200	500	VC
VIS	350-750	145-100*	1800	500	VD
VIS	350-640	75-50*	2400	VIS	VE
NIR	500-1050	500	600	750	NB
NIR	500-1050	220-150*	1200	750	NC
NIR	600-1160	350-300	830	800	SI
NIR	600-1160**	560**	300	1000	IA
NIR	600-1160	500	600	1000	IB

* depends on the starting wavelength of the grating; the higher the wavelength, the bigger the dispersion and the smaller the range to select.

** please note that not all 2048 pixels will be used for the useable range

Resolution Table (FWHM in nm) for AvaSpec-ULS2048x64-EVO

Grating (lines/mm)	Slit size (µm)					
	10	25	50	100	200	500
300	1.40	1.50	2.5	4.8	9.2	21.3
600	0.70 - 0.80*	0.75-0.85*	1.2	2.4	4.6	10.8
830	0.42 - 0.48*	0.50-0.58*	0.93	1.7	3.4	8.5
1200	0.25 - 0.31*	0.37 - 0.43*	0.52-0.66*	1.1	2.3	5.4
1800	0.17 - 0.21*	0.26 - 0.32*	0.34-0.42*	0.8	1.6	3.6
2400	0.12 - 0.18*	0.18 - 0.24*	0.26-0.34*	0.44-0.64*	1.1	2.7
3600	0.09 - 0.12*	0.11 - 0.15*	0.19	0.4	0.8	1.8

* depends on the starting wavelength of the grating; the higher the wavelength, the bigger the dispersion and the better the resolution

Ordering Information

AvaSpec-ULS2048x64-EVO

- Ultra-low Stray-light Fiber-optic Spectrometer, 75 mm AvaBench, 2048x64 pixel back-thinned CCD detector, USB powered, high-speed USB3.0 and ETH interface, incl. AvaSoft-Basic, USB interface cable. Specify grating, wavelength range and options

Options

-RS	• Replaceable slit
DCL-UV/VIS-200	• Quartz detector collection lens (200 - 1100 nm)
SLIT-XX	• Slit size, please specify XX = 5, 10, 25, 50, 100, 200 or 500 µm
SLIT-XX-RS	• Replaceable slit with SMA connector, specify slit size XX = 25, 50, 100, 200 or 500 µm. Only available for AvaSpec-ULS2048CL-EVO-RS
SLIT-XX-RS-FCPC	• As SLIT-XX-RS, but with FC/PC connector
OSF-YYY	• Order-sorting filter for reduction of second-order effects please specify YYY = 305, 395, 475, 515, 550 or 600 nm
OSC	• Order-sorting coating with 600 nm long-pass filter for BB (>305 nm) and VB gratings, recommended with OSF-305
OSC-UA	• Order-sorting coating with 350 and 600 nm linear variable filter for UA, VA gratings
OSC-UB	• Order-sorting coating with 350 and 600 nm long-pass filter for UB or BB (<350 nm) gratings
-FCPC	• FC/PC fiber-optic connector

nordics@amstechnologies.com
azpect.amstechnologies.com

Contact us

Distributor
AZPECT
part of
amstechnologies
where technologies meet solutions

NEW

AvaSpec-ULS2048x64TEC-EVO SensLine Thermoelectrically Cooled Fiber-Optic Spectrometer

The AvaSpec-ULS2048x64TEC-EVO is an updated version of our AvaSpec-ULS2048x64TEC spectrometer, with improved electronics and cooling.

This instrument enhances the Sensline series with its cooled, back-thinned detector. The back-thinned detector has good sensitivity in the UV and IR region. The 64 pixelheight (0.89 mm) enables catching as many photons as possible while the cooling enables long integration times up to 120 seconds with low-noise levels.

The instrument features Peltier cooling device integrated into our exclusive ultra-low stray light optical bench, which can reduce the temperature of the CCD chip to -30°C against ambient, improving the dark baseline and PRNU level significantly. The detector cooling also reduces the dark noise by a factor of 2-3.

The AvaSpec-ULS2048x64TEC-EVO uses a special low-noise version of the 2048x64 detector with integrated cooling.

All the features mentioned above make this instrument ideally suited for measuring low-light applications, such as fluorescence or low-light Raman measurements.

Optimal flexibility is guaranteed with the replaceable slit, making the instrument suitable for various kinds of applications.

The above mentioned qualities make the AvaSpec-ULS2048x64TEC-EVO an excellent choice for low light-level applications, such as fluorescence and Raman measurements, where integration times of more than 5 seconds may be needed.

AvaSpec-ULS2048x64TEC-EVO



Distributor



AZPECT
part of
amSTECHNOLOGIES
where technologies meet solutions

nordics@amstechnologies.com
azpect.amstechnologies.com

Contact us 

Technical Data

Optical bench	ULS Symmetrical Czerny-Turner, 75 mm focal length
Wavelength range	200-1160 nm
Resolution	0.09 –20 nm, depending on configuration (see table)
Stray light	<1%, depending on the grating
Sensitivity	300,000 counts/μW per ms integration time
Detector	Backthinned CCD, 2048x64 pixels, low noise, integrated cooling
Temperature-cooled CCD	Max. ΔT = -30°C versus ambient. Optimal setting: 5°C
Signal/noise	550:1
AD converter	16-bit, 500 KHz
Dynamic range	19,000
Dark noise	5 cnts
Integration time	9.7 ms–120 s
Interface	USB 3.0 high speed, 5 Gbps Gigabit Ethernet 1 Gbps
Sample speed with on-board averaging	9.7 ms/scan
Data transfer speed	9.7 ms/scan (USB3) 9.7 ms/scan (ETH)
Digital IO	HD-26 connector, 2 Analog in, 2 Analog out, 13 Digital bidirectional, trigger, sync., strobe, laser
Power supply	12 VDC, 1.5 A
Operating temperature	0-40°C
Cooling	30°C versus ambient
Dimensions, weight	185 x 145 x 185 mm, 3500 grams

Grating Selection Table for AvaSpec-ULS2048x64TEC-EVO

Use	Useable range (nm)	Spectral range (nm)	Lines/mm	Blaze (nm)	Order code
UV/VIS/NIR	200-1160**	960**	300	300	UA
UV/VIS/NIR	200-1100**	900**	300	300/1000	UNA-DB
UV/VIS	200-850	520	600	300	UB
UV	200-750	250-220*	1200	250	UC
UV	200-650	165-145*	1800	UV	UD
UV	200-580	115-70*	2400	UV	UE
UV	200-400	70-45*	3600	UV	UF
UV/VIS	250-850	520	600	400	BB
VIS/NIR	300-1160**	860**	300	500	VA
VIS	360-1000	500	600	500	VB
VIS	300-800	250-200*	1200	500	VC
VIS	350-750	145-90*	1800	500	VD
VIS	350-640	75-50*	2400	VIS	VE
NIR	500-1050	500	600	750	NB
NIR	500-1050	220-150*	1200	750	NC
NIR	600-1160	350-300	830	800	SI
NIR	600-1160**	560**	300	1000	IA
NIR	600-1160	500	600	1000	IB

* depends on the starting wavelength of the grating; the higher the wavelength, the bigger the dispersion and the smaller the range to select.

** please note that not all 2048 pixels will be used for the useable range

Resolution Table (FWHM in nm) for AvaSpec-ULS2048x64TEC

Grating (lines/mm)	Slit size (µm)					
	10	25	50	100	200	500
300	1.40	1.50	2.5	4.8	9.2	21.3
600	0.70-0.80*	0.75-0.85*	1.2	2.4	4.6	10.8
830	0.42-0.48*	0.50-0.58*	0.93	1.7	3.4	8.5
1200	0.25-0.31*	0.37-0.43*	0.52-0.66*	1.1	2.3	5.4
1800	0.17-0.21*	0.26-0.32*	0.34-0.42*	0.8	1.6	3.6
2400	0.12-0.18*	0.18-0.24*	0.26-0.34*	0.44-0.64*	1.1	2.7
3600	0.09-0.12*	0.11-0.15*	0.19	0.4	0.8	1.8

* depends on the starting wavelength of the grating; the higher the wavelength, the bigger the dispersion and the better the resolution

Ordering Information

AvaSpec-ULS2048x64TEC-EVO

Thermoelectrically cooled fiber-optic spectrometer, 75 mm ultra-low stray light AvaBench, 2048x64 pixel, TE-cooled and regulated low-noise CCD detector, USB3/ETH high-speed interface and replaceable slit, incl. AvaSoft-Basic, USB cable, desktop housing. Specify grating, wavelength range and options

Options

DCL-UV/VIS-200	• Detector Collection Lens to enhance sensitivity, Quartz, 200-1100 nm
SLIT-XX-RS	• Replaceable slit with SMA connector. Specify slit size XX= 10, 25, 50, 100, 200 or 500 µm
SLIT-XX-RS-FCPC	• As SLIT-XX-RS, but with FC/PC connector
OSF-YYY	• Order-sorting filter for reduction of 2 nd order effects, 1 mm thick, please specify YYY= 305, 395, 475, 515, 550 or 600 nm
OSC	• Order-sorting coating with 600 nm long-pass filter for BB (>350 nm) and VB gratings, recommended with OSF-305
OSC-UA	• Order-sorting coating with linear variable filter for UA, VA gratings
OSC-UB	• Order-sorting coating with 350 and 600 nm long-pass filter for UB or BB (<350 nm) gratings

nordics@amstechnologies.com
aspect.amstechnologies.com

Contact us

AvaSpec-ULS2048XL-EVO SensLine High UV and NIR Sensitivity Back-thinned CCD Spectrometer

AvaSpec-ULS2048XL-EVO



Combining exceptional quantum efficiency with high-speed is the value proposition of the AvaSpec-ULS2048XL-EVO spectrometer. Unlike many back-thinned CCD spectrometers, which have two dimensional arrays, the ULS2048XL-EVO has large monolithic pixels of 14x500 microns with exceptional efficiency in the UV, from 200-400 nm, and the NIR, from 950-1160 nm. The instrument also has an electronic shutter, which enables integration times as low as 2 microseconds. To further enhance sensitivity, a detector collection lens is available which improves sensitivity up to 60% when combined with larger core fibers.

Options include order-sorting filter, to reduce 2nd order effects and purge ports for deep-UV measurements. The AvaSpec-ULS2048XL-EVO comes with a wide range of slit sizes, gratings and may be configured with SMA or FC/PC fiber-optic entrance connectors.

The AvaSpec-ULS2048XL-EVO uses the AS7010 electronics board offering USB3 (10 times faster than USB2), Gigabit Ethernet and better signal processing.

Connection to your PC is handled via a USB3-connection or Ethernet, delivering a scan every 2 milliseconds. The instrument comes complete with AvaSoft-basic software, USB cable and an extensive manual.

Technical Data

Optical Bench	ULS, Symmetrical Czerny-Turner, 75 mm focal length
Wavelength range	200 - 1160 nm
Resolution	0.09 -20 nm, depending on configuration (see table)
Stray-light	< 0.5%
Sensitivity	460,000 counts/ μ W per ms int. time
UV Quantum efficiency	60% (200-300 nm)
Detector	Back-thinned CCD image sensor 2048 pixels
Signal/Noise	525:1
AD converter	16-bit, 1 MHz
Integration time	2 μ s - 20 seconds
Interface	USB 3.0 high-speed, 5 Gbps Gigabit Ethernet 1 Gbps
Sample speed with store to RAM	2.44 ms /scan
Readout Noise	9.8 cnt RMS
Dark Noise	4.5 cnt RMS
Dynamic Range	13.700
Data transfer speed	2.44 ms /scan (USB3)
Digital IO	HD-26 connector, 2 Analog in, 2 Analog out, 3 Digital in, 12 Digital out, trigger, synchronization
Power supply	Default USB power, 700 mA. Or external 12VDC, 360 mA
Dimensions, weight	175 x 127 x 44,5 mm (1 channel), 1180 grams

Distributor



AZPECT
part of
amSTECHNOLOGIES

nordics@amstechnologies.com
aspect.amstechnologies.com

where technologies meet solutions

Contact us 

AVANTES
MEMBER OF THE NYNOMIC GROUP

Grating Selection Table for AvaSpec-ULS2048XL-EVO

Use	Useable range (nm)	Spectral range (nm)	Lines/mm	Blaze (nm)	Order code
UV/VIS/NIR	200-1160**	960**	300	300	UA
UV/VIS/NIR	200-1160**	900**	300	300/1000	UNA-DB
UV/VIS	200-850	520	600	300	UB
UV	200-750	250-220*	1200	250	UC
UV	200-650	165-145*	1800	UV	UD
UV	200-580	115-70*	2400	UV	UE
UV	200-400	70-45*	3600	UV	UF
UV/VIS	250-850	520	600	400	BB
VIS/NIR	300-1160**	860**	300	500	VA
VIS	360-1000	500	600	500	VB
VIS	300-800	250-200*	1200	500	VC
VIS	350-750	145-100*	1800	500	VD
VIS	350-640	75-50*	2400	VIS	VE
NIR	500-1050	500	600	750	NB
NIR	500-1050	220-150*	1200	750	NC
NIR	600-1160	350-300	830	800	SI
NIR	600-1160**	560**	300	1000	IA
NIR	600-1160	500	600	1000	IB

* depends on the starting wavelength of the grating; the higher the wavelength, the bigger the dispersion and the smaller the range to select.

** please note that not all 2048 pixels will be used for the useable range

Resolution Table (FWHM in nm) for AvaSpec-ULS2048XL-EVO

Grating (lines/mm)	Slit size (μm)					
	10	25	50	100	200	500
300	1.40	1.50	2.5	4.8	9.2	21.3
600	0.70 - 0.80*	0.75-0.85*	1.2	2.4	4.6	10.8
830	0.42 - 0.48*	0.50-0.58*	0.93	1.7	3.4	8.5
1200	0.25 - 0.31*	0.37 - 0.43*	0.52-0.66*	1.1	2.3	5.4
1800	0.17 - 0.21*	0.26 - 0.32*	0.34-0.42*	0.8	1.6	3.6
2400	0.12 - 0.18*	0.18 - 0.24*	0.26-0.34*	0.44-0.64*	1.1	2.7
3600	0.09 - 0.12*	0.11 - 0.15*	0.19	0.4	0.8	1.8

* depends on the starting wavelength of the grating; the higher the wavelength, the bigger the dispersion and the better the resolution

Ordering Information

AvaSpec-ULS2048XL-EVO

PS-12V / 1.0A

- Ultra-low Stray-light Fiber-optic Spectrometer, 75 mm AvaBench, 2048 large 500 μm pixel back-thinned CCD detector, USB powered, high-speed USB3.0 and ETH interface, incl. AvaSoft-Basic, USB interface cable.
Specify grating, wavelength range and options
- External power supply, needed for use in ETH mode

Why is the XL so sensitive?
We're using back-illuminated detectors.
They have the electronics on the backside of the detector,
allowing more light to be caught by the front side.

Options

-RS	• Replaceable slit
DCL-UV/VIS-200	• Quartz Detector Collection Lens (200-1100 nm)
SLIT-XX	• Slit size, please specify XX = 10, 25, 50, 100, 200 or 500 μm
SLIT-XX-RS	• Replaceable slit with SMA connector , specify slit size XX=25, 50, 100, 200 or 500 μm . Only in combination with AvaSpec-ULS2048XL-EVO-RS
SLIT-XX-RS-FCPC	• as SLIT-XX-RS, but with FC/PC connector
OSF-YYY	• Order-sorting filter for reduction of 2nd order effects, 1 mm thick, please specify YYY= 305, 395, 475, 515, 550 or 600 nm
OSC	• Order-sorting coating with 600 nm long-pass filter for BB (>350 nm) and VB gratings, recommended with OSF-305
OSC-UA	• Order-sorting coating Linear Variable Filter for UA, VA gratings
OSC-UB	• Order-sorting coating with 350 and 600 nm long-pass filter for UB or BB (<350 nm) gratings
-FCPC	• FC/PC fiber-optic connector

Distributor



amstechnologies
where technologies meet solutions



AZPECT
part of

nordics@amstechnologies.com
azpect.amstechnologies.com

Contact us 

The grating can only be changed by Avantes.
Therefore, choose your grating wisely.
Our application specialists are available to support you with your choice.
In general, a higher resolution means a lower bandwidth.
By combining multiple spectrometers
in our AvaSpec-Dual or rack-mountable versions,
you can create one virtual spectrometer with high-resolution
and high bandwidth.